

# MATERIAL SAFETY DATA SHEET



Date Issued: 08/31/2010  
MSDS No: 2010.263

## 1. PRODUCT AND COMPANY IDENTIFICATION

**PRODUCT CODE:** E6380  
**PRODUCT FORMULATION NAME:** Ethyltriethoxysilane  
**MOLECULAR FORMULA:** C<sub>8</sub>H<sub>20</sub>O<sub>3</sub>Si

### MANUFACTURER

UCT Specialties, Inc.  
2731 Bartram Road  
Bristol, PA 19007

**Emergency Contact:** Jon Telepchak  
**E-Mail:** jtelepchak@unitedchem.com  
**Product Stewardship:** 215-781-9255  
ex141  
**Service Number:** 717-247-0896

### 24 HR. EMERGENCY TELEPHONE NUMBERS

**CANUTEC (Canadian Transportation) :**(613) 996 - 6666  
**CHEMTREC (US Transportation) :**(800) 424 - 9300

## 2. HAZARDS IDENTIFICATION

### EMERGENCY OVERVIEW

**PHYSICAL APPEARANCE:** Clear liquid

**IMMEDIATE CONCERNS:** Material may form a siloxane polymer on the skin, eyes or in the lungs. In the event of direct contact of the liquid with these tissues, seek medical attention.

### POTENTIAL HEALTH EFFECTS

**EYES:** Severely irritating. Conjunctivitis, corneal damage.

**SKIN:** \*Irritant. May be harmful. Avoid contact .

**INGESTION:** \*Irritant. Effects undetermined. Treat as toxic

**INHALATION:** \* Causes irritation or damage to lung

### SIGNS AND SYMPTOMS OF OVEREXPOSURE

**ACUTE TOXICITY:** Material generates Ethanol on contact with water or moisture on skin, eyes, and mucous membranes and has an irritating, dehydrating effect on overexposed tissue.

Overexposure to Ethanol by skin absorption, inhalation, or ingestion may have a narcotic effect (headache, nausea, drowsiness). Ethanol is metabolized to acetaldehyde and acetic acid resulting in metabolic acidosis, CNS depression, and death due to respiratory arrest.

## 3. COMPOSITION / INFORMATION ON INGREDIENTS

Chemical Name	CAS
Ethyltriethoxysilane	78-07-9

## 4. FIRST AID MEASURES

**EYES:** Flush with clean water for at least 15 minutes and consult physician

**SKIN:** Scrub with soap and water. Remove contaminated clothing and shoes

**INGESTION:** Give one full cup of milk or water to dilute ingested material. Get medical attention

**INHALATION:** Remove victim to fresh air, give CPR or oxygen if necessary

**COMMENTS:** If symptoms persist, get medical attention

Never give anything by mouth to an unconscious person

## 5. FIRE FIGHTING MEASURES

**FLASH POINT AND METHOD:** 40°C open cup

**EXTINGUISHING MEDIA:** Water, CO<sub>2</sub>, foam, dry chemical

**FIRE FIGHTING EQUIPMENT:** Fire Fighters must wear positive- pressured, self-contained breathing apparatus and full protective clothing.

## 6. ACCIDENTAL RELEASE MEASURES

**GENERAL PROCEDURES:** Absorb spilled material with suitable chemical binder. Shovel absorbent into suitable waste container. Flush spill area with large amount of water. Do not contaminate soil, groundwater, or surface water.

**SPECIAL PROTECTIVE EQUIPMENT:** Wear appropriate personal protective equipment as specified in section 8.

## 7. HANDLING AND STORAGE

**HANDLING:** Container requires grounding during use.

**STORAGE:** Store in a cool and dry place. Protect from moisture. Protect from heat, direct sunlight, and source of ignition. Maintain tightly closed container. Store away from alkaline, acidic, and oxidizing materials. Provide adequate ventilation

## 8. EXPOSURE CONTROLS / PERSONAL PROTECTION

**ENGINEERING CONTROLS:** Facilities storing or utilizing this material should be equipped with an eyewash facility and a safety shower. Local exhaust ventilation may be necessary to control any fume levels during the use of this product.

### PERSONAL PROTECTIVE EQUIPMENT

**EYES AND FACE:** Use Chemical goggles.

**SKIN:** Use impervious gloves. Use impervious clothing as necessary to protect against skin contact.

**RESPIRATORY:** Maintain adequate ventilation, Mechanical ventilation recommended

In case of exposure, use appropriate NIOSH approved respiratory protection

**WORK HYGIENIC PRACTICES:** General industrial hygiene practice.

**COMMENTS:** To identify additional Personal Protective Equipment (PPE) requirements, it is recommended that a hazard assessment in accordance with the OSHA PPE Standard (29 CFR 1910.132) be conducted before using this product.

## 9. PHYSICAL AND CHEMICAL PROPERTIES

**PHYSICAL STATE:** Liquid

**COLOR:** Clear

**VAPOR DENSITY:** > 1

**BOILING POINT:** 158°C /760mm

**FLASH POINT AND METHOD:** 40°C open cup

**SOLUBILITY IN WATER:** Reacts

**SPECIFIC GRAVITY:** 0.896

## 10. STABILITY AND REACTIVITY

**STABILITY:** (UNDER NORMAL CONDITIONS): Stable under ordinary conditions of use and storage.

**POLYMERIZATION:** Hazardous polymerization will not occur

**CONDITIONS TO AVOID:** Contact with water and moist air ; Generates Hydrogen Chloride;  
TVL=C5ppm. Exposure to heat

**INCOMPATIBLE MATERIALS:** oxidizing agents, alkali and acid

## 11. TOXICOLOGICAL INFORMATION

### ACUTE

**SKIN ABSORPTION:** 16 g/kg- rbt

**ORAL LD<sub>50</sub>:** 14 g/kg-rat

### CARCINOGENICITY

**IARC:** No

**NTP:** No

**OSHA:** No

## 12. ECOLOGICAL INFORMATION

**COMMENTS:** No information available

## 13. DISPOSAL CONSIDERATIONS

**DISPOSAL METHOD:** Incineration is recommended . Treat or dispose of waste material in accordance with all local, state/provincial and national requirements.

**EMPTY CONTAINER:** Containers of this material may be hazardous when empty since they retain product residues(liquid); observe all warnings and precautions listed for this product.

## 14. TRANSPORT INFORMATION

### DOT (DEPARTMENT OF TRANSPORTATION)

**PROPER SHIPPING NAME:** Flammable liquids, n.o.s.

**TECHNICAL NAME:** (Triethxy(ethyl)silane)

**PRIMARY HAZARD CLASS/DIVISION:** 3

**UN/NA NUMBER:** UN1993

**PACKING GROUP:** III

## 15. REGULATORY INFORMATION

### UNITED STATES

#### DOT LABEL SYMBOL AND HAZARD CLASSIFICATION



Flammable Liquid

### TSCA (TOXIC SUBSTANCE CONTROL ACT)

**TSCA STATUS:** This product is listed on the TSCA Inventory.

## 16. OTHER INFORMATION

**APPROVED BY:** Office of Environmental Health Safety & Security  
Security

**TITLE:** Dir. of Env. Health Safety &

**PREPARED BY:** jnm


**INFORMATION CONTACT:** jtelepchak@unitedchem.com

**REVISION SUMMARY:** New MSDS

#### HMIS RATING

<b>HEALTH:</b>	<b>2</b>
<b>FLAMMABILITY:</b>	<b>2</b>
<b>PHYSICAL HAZARD:</b>	<b>1</b>
<b>PERSONAL PROTECTION:</b>	

#### NFPA CODES

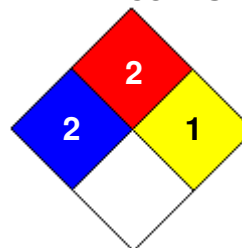

**MANUFACTURER SUPPLEMENTAL NOTES:** All technical data is believed to be correct but does not purport to be all inclusive and shall be used only as a guide.

**MANUFACTURER DISCLAIMER:** The data in this Material Safety Data Sheet relates only to the specific material designated herein.

It does not relate to use in combination with any other material or in any process.

THE INFORMATION AND RECOMMENDATIONS CONTAINED IN THIS MATERIAL SAFETY DATA SHEET HAVE BEEN COMPILED FROM SOURCES BELIEVED TO BE RELIABLE AND TO REPRESENT THE BEST OPINION ON THE SUBJECT AS OF THE DATE ON THIS SHEET. HOWEVER, NO WARRANTY, GUARANTEE OR REPRESENTATION, EXPRESS OR IMPLIED, IS MADE BY UCT, INC. AS TO THE CORRECTNESS OR SUFFICIENCY OF THIS INFORMATION OR TO THE RESULTS TO BE OBTAINED FROM THE USE THEREOF.