




QUICKLY PROCESS UP TO 96 SAMPLES AT ONE TIME

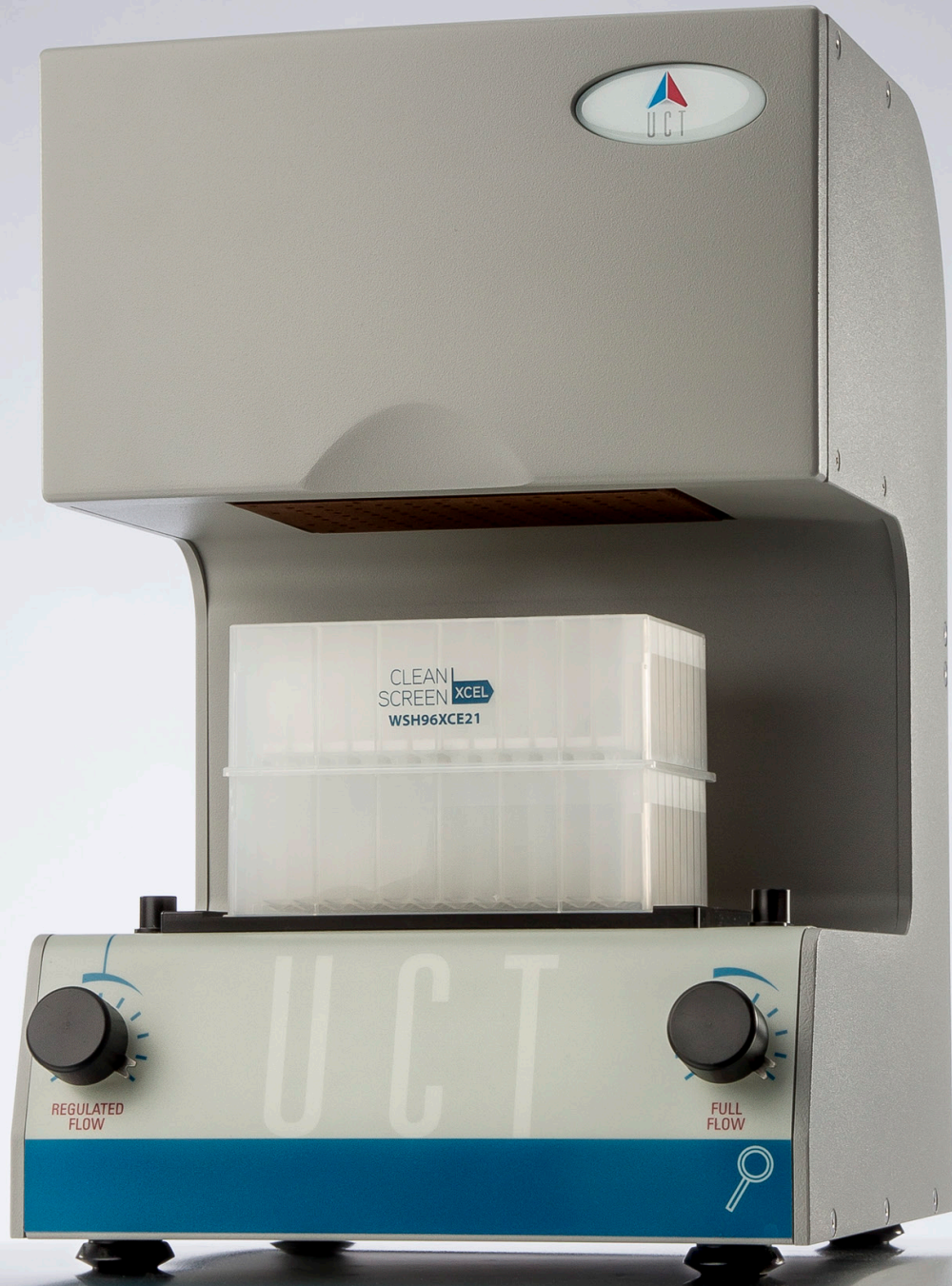

96 WELL PLATE POSITIVE PRESSURE MANIFOLD

96 WELL PLATE POSITIVE PRESSURE MANIFOLD

- Small instrument footprint (9" x 9" x 14") does not take up valuable bench space
- Accommodates a wide range of commercially available 96 well plates
- Works with both 96 and 48 well plates
- Restrictor plate allows even gas flow across all 96 wells – even when some positions are empty
- Two gas flow regulators on the front panel – one high flow regulator used for sorbent bed drying and one low flow regulator used for sample loading, washing, and eluting
- Positive pressure is well suited for viscous sample matrices
- Instrument gas regulator on the back can be set and locked into place
- Gas Supply: N₂ or compressed air regulated to 40 – 45 psi and filtered to 10µm
- All pneumatic – only a gas supply is needed, there are no electrical connections

PPM		
Part Number	Description	Unit
VMF96PPM	96 Well Plate Positive Pressure Manifold	Ea
*For use of complete Positive Pressure Manifold System, PPM installation kit must also be purchased.		
Accessories		
Part Number	Description	Unit
VMFPPMIK	Installation kit (25' of ¼" O.D. tubing; 2 x ¼" compression fittings; 1 in-line air filter)	Kit
VMF96PPMGSK	Replacement Brown Gasket – 96 well plate manifold	Ea
WSH96WT	96 well waste collection plate	Ea
WSH96CP	96 well sample collection plate	Ea

Contact Us

Phone: 215.781.9255
800.385.3153
Fax: 215.785.1226

UCT, Inc.
2731 Bartram Rd.
Bristol, PA 19007

Email: info@unitedchem.com
Web: www.unitedchem.com



4101-18-04

