


Sympathomimetic Amines In Blood, Plasma/Serum, Urine, and Tissue For GC Or GC/MS Confirmations Using: 200 mg Clean Screen® Extraction Column

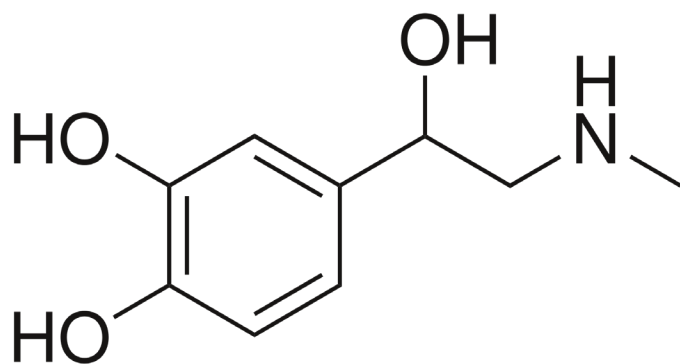

UCT Part Numbers

ZSDAU020
Clean Screen® DAU
10 mL, 200 mg Column

Or

ZCDAU020
Clean Screen® DAU
10 mL, 200 mg Column
with CLEAN-THRU® Tips

Procedure:

1. Prepare Sample

- To 1 mL of 100 mM phosphate buffer (pH 6.0) add internal standard(s).^{*} Add 1 mL of blood, plasma/serum, urine or 1g of (1:4) tissue homogenate.
- Mix/vortex.
- Add 2 mL of 100 mM phosphate buffer (pH 6.0). Mix/vortex.
- Sample pH should be 6.0 ± 0.5 .
- Adjust pH accordingly with 100 mM monobasic or dibasic sodium phosphate.
- Centrifuge as appropriate.

2. Condition Clean Screen® Extraction Column

- 1 x 3 mL CH₃OH.
- 1 x 3 mL D.I. H₂O.
- 1 x 1 mL 100 mM phosphate buffer (pH= 6.0).

Note: Aspirate at < 3 inches Hg to prevent sorbent drying.

3. Apply Sample

- Load at 1 mL/minute.

4. Wash Column

- 1 x 3 mL D.I. H₂O.
- 1 x 1 mL 100 mM acetic acid.
- 1 x 3 mL CH₃OH.
- Dry column (5 minutes at > 10 inches Hg).

5. Elute SMA

- 1 x 3 mL CH₂Cl₂/IPA/NH₄OH (78:20:2); Collect eluate at 1 to 2 mL/minutes.

Note: Prepare elution solvent daily. Add IPA/NH₄OH, mix, then add CH₂Cl₂

6. Concentrate Eluate

- Add 30 µL silylation grade DMF to eluate.
- Evaporate to 30 µL at < 40°C.

Alternate Drying Procedure

- Evaporate for 4 min.
- Add 100 µL of 1% HCl in methanol.
- Evaporate to dryness.

7. Fluoroacylate With PFPA (PFAA)***

- Add 50 µL PFPA (PFAA). Overlay with N² and cap.
- Improved derivatization by addition of 50 µL PFPOH.****
- React 20 minutes at 70°C. Evaporate to dryness at < 40°C.
- Reconstitute with 100 µL ethyl acetate.

8. Quantitate

- Inject 1 to 2 µL onto gas chromatograph.
- For MSD monitor the following ions: (Next Page)



Analyte (PFPA)	Primary Ion **	Secondary	Tertiary	Cerilliant #
Amphetamine-D5*	194	92	123	A-005
Amphetamine	190	91	118	A-007
Methamphetamine-D5*	208	92	163	M-004
Methamphetamine	204	91	160	M-009
Pseudoephedrine	204	160	119	P-035
Ephedrine	204	160	119	E-011
Phenylephrine	190	119	267	
Methylenedioxyamphetamine	135	162	325	M-012
Methylenedioxymethamphetamine	204	162	339	M-013

* Suggested internal standards for GC/MS: D₅-Amphetamine and D₅-Methamphetamine

** Quantification Ion

*** Part # SPFAA-0-1,10, 25,100

**** Part # SPFPOH-1,10,25,100

Alternate Derivatization

7. Form TMS derivatives:

- Add 50 µL BSTFA with 1% TMCS and 50 µL of ethyl acetate.
- React 45 minutes at 70°C.

8. Quantitate

- Inject 1 to 2 µL onto gas chromatograph.
- For MSD monitor the following ions:

Analyte (TMS)	Primary Ion **	Secondary	Tertiary	Cerilliant #
Amphetamine-D5*	120	197	92	A-005
Amphetamine-D6*	120	198	93	A-044
Amphetamine-D10*	120	202	97	A-038
Amphetamine-D12*	120	203	98	A-019
Amphetamine	116	192	91	A-007
Methamphetamine-D5*	134	211	92	M-023
Methamphetamine-D8*	137	214	92	
Methamphetamine-D9*	137	215	93	M-091
Methamphetamine	130	206	91	M-004
Pseudoephedrine	130	147	294	P-035
Ephedrine	130	147	294	E-011
Methylenedioxyamphetamine	116	236	135	M-012
Methylenedioxymethamphetamine	130	250	131	M-013
Para-Methoxamphetamine	116	222	121	NMID1908



Alternate Derivatization

7. Form 4-CB:

- (4-Carbethoxyhexafluorobutyl chloride)* derivatives:
Add 20 µL 4-CB* and 100 µL of ethyl acetate.
- React 45 minutes at 70°C.

8. Quantitate

- Inject 1 to 2 µL onto gas chromatograph.
- For MSD monitor the following ions:

Analyte (TMS)	Primary Ion **	Secondary	Tertiary	Cerilliant #
Amphetamine-D5****	298	270	399	A-005
Amphetamine	294	266	248	A-007
Methamphetamine-D5****	312	284	266	M-023
Methamphetamine	308	280	262	M-004
Methylenedioxyamphetamine-D5****	136	434	270	M-010
Methylenedioxyamphetamine	162	429	266	M-012
Methylenedioxymethamphetamine-D5****	312	284	266	M-011
Methylenedioxymethamphetamine	308	280	262	M-013
Methylenedioxyethylamphetamine-D6****	328	165	300	
Methylenedioxyethylamphetamine	322	162	294	M-065

* Part # S4CB-0-10

*** Quantification Ion

**** Suggested internal standards for GC/MS: D5-Amphetamine and D5-Methamphetamine

UCT, LLC • 2731 Bartram Road • Bristol, PA 19007 800.385.3153 • 215.781.9255

www.unitedchem.com Email: methods@unitedchem.com

©UCT, LLC 2023 • All rights reserved

